

Evolve Hypromellose

0.3%

Hypromellose 0.3% Eye Drops
Preservative free 10ml

Evolve Hypromellose is a sterile ophthalmic solution for adults (18-64 years) with physiological pH, based upon hydroxypropyl methylcellulose.

The main action of the device is to relieve dry eye symptoms by moisturising and lubricating the ocular surface thanks to hydroxypropyl methylcellulose. It soothes and reduces dry eye symptoms such as burning, itching, irritation and foreign body sensation, and protects the ocular surface by stabilising the tear film and creating a protective, transparent and viscoelastic shield.

Evolve Hypromellose relieves eye dryness due to digital device use, infrequent blinking, environmental factors and wearing contact lenses.

Evolve Hypromellose is compatible with contact lenses.

INTENDED PURPOSE

To relieve eye discomfort, burning, itching and irritation due to environmental factors of different origins.


INDICATIONS

It is indicated in case of moderate dry eye syndrome, dry eye symptoms or discomfort sensation.

INSTRUCTIONS FOR USE

- Wash your hands carefully before instilling the product.
- At first use, remove the tamper-evident ring that seals the cap by pulling off the side tab (figure 1).
- Remove the cap (figure 2).
- Place the cap taking care to put it on a clean surface.
- Turn the bottle upside down over the eye, gripping it in the middle with your thumb and index finger (figure 3).
- Apply an increasing pressure on the sides until the drop comes out (figure 4).
- Instill 1-2 drops in each eye. Wait a few seconds after each drop without squeezing the bottle. This will enable enough air to flow back into the system, otherwise it is more difficult to get the drops out.

- Shake off the bottle to remove any residual drop (figure 5) and then the cap can be put back on (figure 6). If, after use, the bottle is squeezed, wait for it to return to its original shape before closing it with the cap.

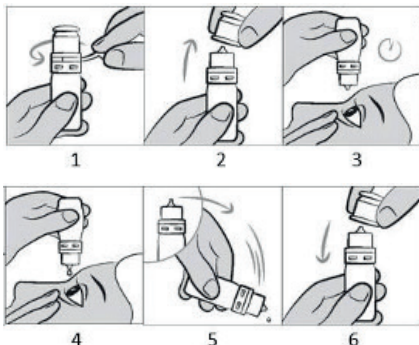

- The product is also indicated for frequent use, even several times a day (3 or 4 times) or according to healthcare professional indication.

COMPOSITION

Hydroxypropyl methylcellulose (HPMC) 0.3%, Sodium phosphate monobasic monohydrate, Sodium phosphate dibasic dodecahydrate, Sodium chloride, Water for injections.
It does not contain preservatives.

WARNINGS AND PRECAUTIONS

- The content must be used within 3 months after the first opening.
- If you suffer from severe damage to the clear layer at the front of the eye (the cornea), phosphates may cause, in very rare cases, cloudy patches on the cornea due to calcium build-up during treatment.
- Follow the prescribed posology or recommended use.
- For topical ophthalmic use only.
- The product is for single patient use.
- At the time of first opening, do not use the product if container is damaged or not perfectly closed.
- Do not use if there are visible particles and/or changes in the colour of the solution.
- Avoid contact between the tip of the bottle and your eyes, fingers or any other surface.
- Do not damage or manipulate or disassemble the container e.g. by the penetration of the orifice, mechanical impact or misuse.

- Do not damage or block the orifice and ventilation hole.
- Do not refill or reuse the pack.
- Do not squeeze the device with attached cap.
- Actuate the device only in tip down position.
- Remove the cap for application and put it back on immediately after using the product.
- The gradual emptying of the bottle may be associated with the need to exert higher pressure to allow the delivery of the product. This is due to the properties of the closure system act to preserve the sterility and the chemical and physical characteristics of preservative free solutions. A small quantity of solution also always remains in the bottle.
- Do not use the product after the expiry date and the period after first opening indicated on the package.
- Please ensure that the recommended storage conditions are always observed.
- The product can be used for up to 30 consecutive days. If symptoms persist beyond this period, you should consult your healthcare professional.
- In case of side effects during treatment, stop using this product and consult a healthcare professional.
- Do not use in case of known incompatibility towards any of the components of the solution.
- No interactions with drugs or other eye drops are reported; however, concomitant use with ophthalmic drugs should be avoided and at least 15 minutes should be allowed between administration of the product and any other eye drops.
- After use, dispose of the product according to your local waste management regulations.
- Keep out of reach and sight of children. It is recommended to consult your healthcare professional in case of use on children.
- Pregnant or breast-feeding women and those who have recently undergone eye surgery must consult their healthcare professional before use.
- This product does not substitute for any drug that may be prescribed by your healthcare professional.
- For external use, do not swallow.
- It is necessary to report any serious incident occurred in relation to the device to the manufacturer and the competent authority of the member state in which you are established.
- The product is not indicated for the treatment of eye infections and other ocular diseases. These conditions require medical attention.

SIDE EFFECTS

Vision may be temporarily blurred after use of product. Also, minor burning/stinging/irritation may temporarily occur as a consequence of the ocular irritation. If any of these effects persist or worsen, tell your healthcare professional promptly.

CONTENTS OF THE PACK

One bottle containing 10 ml of ophthalmic solution.

STORAGE CONDITIONS

Store below 30°C.

VALIDITY

With unopened pack: do not exceed the expiry date stated on the carton and bottle.

Last revision of the text: 02.01 of 03/2026

The Evolve brand and logo are trademarks of Medicom Healthcare Limited.



DT/01/DT/015

CE 0051

 C.O.C. Farmaceutici S.r.l.
Via Modena, 15
40019 - Sant'Agata Bolognese (BO) Italy

 Medicom Healthcare Limited,
Lynton House, 7-12 Tavistock Square,
Kings Cross, London, WC1H 9LT, U.K.