

Evolve® Revive

**Sodium Hyaluronate 0.2% w/v +
Carbomer 980 0.2% w/v**

**Preservative Free
10ml**

Instructions for use

It is important to read this leaflet carefully before using **Evolve® Revive**. The information provided summarises the important points about the product but if you have any further questions then please ask your doctor, pharmacist or optician.

LUBRICATING, PROTECTING
& HYDRATING EYE DROPS

23373802

EASY SQUEEZE BOTTLE

Soft bottle to improve
ease of use



1. What is Evolve® Revive?

Evolve® Revive is a multi-dose, preservative free, sterile solution which contains Sodium Hyaluronate 0.2% w/v + Carbomer 980 0.2% w/v in a 10ml bottle. The Evolve® system has been designed to combine the benefits of unit dose preservative free eye drops with the convenience of a standard eye drop bottle.

Evolve® Revive contains Sodium Hyaluronate, Carbomer 980 and Glycerol. Each has a distinct characteristic to provide relief from dry eyes sensations. Sodium Hyaluronate is a naturally occurring, biocompatible material with unique viscoelastic and hygroscopic properties while the Carbomer 980 provides long-lasting relief and the Glycerol acts as a humectant and protects the eye surface.

The combination provides a Triple Action on the eyes surface.

2. Composition:

Evolve® Revive contains Sodium Hyaluronate 0.2% w/v, Carbomer 980 0.2% w/v, Glycerol, Sorbitol, Sodium Chloride, Water for Injections and a pH modifier.

Evolve® Revive is preservative free, phosphate & buffer free ideal for frequent and occasional users.

It has been demonstrated that all eyes should not be exposed to these irritants.

3. What is Evolve® Revive used for?

Evolve® Revive is a lubricating eye drop intended for use in the relief of discomfort that arises from dry eye sensations. It provides long-lasting relief for gritty, itchy and irritated dry eyes and other foreign body sensations.

The lack of irritants such as buffers and preservatives, make this product suitable for prolonged and protective usage in patients with sensitive eyes including itchy and excessive watering of the eyes associated with an unstable tear film.

Dry eye sensations may also be produced by the use of smartphones, TVs and computer screens, infrequent blinking, certain medical treatments including refractive laser eye surgery, atmospheric pollution and drying atmosphere e.g. air conditioning, central heating, wind and sun.

Evolve® Revive is compatible with contact lenses.

The bottle has been designed to assist ease of use with patients. The body of the bottle is easy to grip using a pinch grip between the thumb and forefinger and has been constructed of a soft pliable plastic making squeezing and control of dispensing a drop easy and convenient. Only GENTLE pressure is required to squeeze one drop out. The blue tip contrasts against a white background, that will assist in aiming the drop as you invert the bottle towards your eye. Lastly, there is no protruding spout that could poke and damage your eye.

The bottle delivers 360 calibrated drops of solution.

4. How to use:

Before opening a new bottle of **Evolve® Revive**, please check the safety seal of the bottle is unbroken by checking the integrity of the six bridges between the tamper-evident ring and the cap.


- Wash and dry your hands well before use.
- Rotate the cap anticlockwise to break the seal.
- Remove the cap.
- Tilt your head back and look up.
- Hold the bottle between your thumb and forefinger of your hand above your eye.
- Gently pull down the lower lid of the eye to form a small pocket between your eyelid and your eye.
- Apply moderate pressure in the middle of the bottle to release one drop of **Evolve® Revive** into



the pocket. Please note that there may be a few seconds delay between squeezing and the drop coming out. **Do not squeeze too hard.**

- Blink a few times to spread the drop over the eye.
- If you think that you have missed the eye, then insert another drop.
- Carefully wipe away any excess product around the eyelid with a tissue.
- Repeat the process for the other eye, if required.
- Shake the bottle once downwards to remove any excess drop (see Figure 1) before replacing the cap on the container, do not touch the applicator tip with anything, including your eye or your fingers.

If you are unsure about how to use **Evolve® Revive**, consult your doctor, pharmacist or optician.

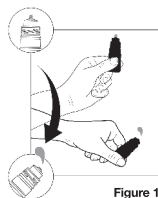

5. How much to use?

The usual dosage is one drop to be applied two times a day as needed, or as required depending on the severity of your dry eye sensation. Your healthcare practitioner may advise you to use this product in a different way, so follow their directions about when and how to use this product and read the label.

6. Precautions and warnings:

- This product is for external ophthalmic use only.
- Do not use if the safety seal on the cap of the bottle is broken before first use.
- Do not use if you are allergic to any of the ingredients.
- Do not touch the tip of the bottle with your eye or any surface or dip the nozzle in any liquid.
- Do not share this product with any other person.
- This product can cause a mild burning sensation and/or blurring of vision after application. Wait until your vision is clear before driving or using machinery.
- Do not use after the expiry date shown on the carton or printed on the bottle.
- Leave at least 5 minutes between administrations of any other ophthalmic drugs.
- Keep out of the reach and sight of children.
- Do not deliver more than 30 consecutive drops over a short period.
- Do not try to unscrew or disassemble the nozzle from the bottle.
- Do not attempt to refill the bottle.
- Do not squeeze the bottle when it is in an upright position.
- Consult your doctor, pharmacist or optician before you use this product if you are pregnant or breast-feeding or have undergone an ophthalmic procedure.
- If your eye condition persists, consult your doctor, pharmacist or optician.

7. Contraindications

Do not use if you are allergic to Carbomer and Sodium Hyaluronate or to any of the ingredients listed in section 2-Composition.

8. Possible Side Effects

- Transient blurring of vision
- Mild discomfort in eye
- possible allergic reaction

If you experience any additional adverse events, consult your healthcare professional and report to Medicom Healthcare Ltd.

9. How to store Evolve® Revive?

Do not store above 25°C.

Discard the **Evolve® Revive** three months after first opening.

Dispose of the bottle according to your local waste management regulations.

Bottle Discards Date Reminder: Patients using this product are advised to fill out their name, first open date & discard date in the designated area on the top inside flap of the carton box, this feature will enable the patient to know exactly when to discard the product after first opening.

This leaflet applies to **Evolve® Revive** and not to any other product.

Last revision of the text: May 2023.



Discard the contents 3 months after first opening



Do not use if package is damaged and consult instructions for use



Sterilised using aseptic processing techniques



Consult instructions for use



Do not store above 25°C. Do not freeze the product



Medical device symbol



Medicom Healthcare Ltd.,
Lynton House,
7-12 Tavistock Square, Kings Cross,
London, WC1H 9LT, U.K.

UK
CA
0120