

PACKAGE INSERT




Hypromellose Eye Drops, 0.3% w/v

Hypromellose 3.2 mg
Water for Injection q.s.

Every Hypromellose Eye Drops unit contains sterile 0.3% w/v Hypromellose free from any preservative, for single use.


Intended use:

The Hypromellose Eye Drops, 0.3% w/v is intended to use for the symptomatic lubrication of “dry eye” which is a dehydration of the surface of the eye because no natural tear are being produced.

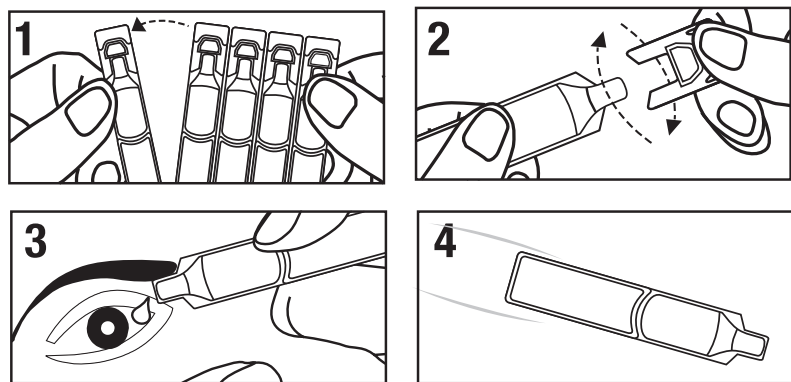

Instructions for use:

1. Detach the container to be used.
2. To open, twist and pull the tab to remove.
3. Instill 1 or 2 drops in the eye(s) and discard the container.
4. Not for multiple use. Discard the vial after single use.



Storage:

- Store between 15°C-25°C
- Do not Freeze
- Protect from (Sun) light.
- Keep out of reach of children

Precautions:

- Do not use, if vial is damaged.
- Do not use, if solution is cloudy or discoloured.
- Do not touch the eye with the vial; this is to avoid eye infection or scratching.
- Not for treatment of any pathological eye conditions.
- Not to be used after the expiry date, marked on the package.
- Dispose of the vial as per the Country specific requirements.
- This product should not be re-used as it may cause microbial infection.

PACKAGE INSERT



Hypromellose Eye Drops, 0.3% w/v

- Do not resterilize the final packed product as it is not suitable for resterilization.

Stop use and ask a doctor if you experience

- Eye pain.
- Changes in vision.
- Continued redness or irritation of the eye, or if the condition worsens or persists for more than 72 hours.

Date of last revision of text: Sept 2025

NOT FOR INJECTION

MD

STERILE A

Sterilised using aseptic processing techniques

15°C – 25°C
Store between 15°C and 25°C



Do not resterilise



Caution, see instructions for use



Do not use if package is damaged



Keep Dry



Do not reuse



Consult instructions for use



Keep away from Sunlight



Sentiss Pharma Pvt. Ltd.
Vill. Khara Nihla, Teh. Nalagarh,
Dist. Solan (HP)-174101, INDIA

CE 2460

EU REP

OBELIS S.A
Boulevard General Wahis 53
1030 Brussels, BELGIUM
Tel: +(32) 2.732.59.54
Fax: +(32) 2.732.60.03
E-Mail: mail@obelis.net



Imported and distributed by:
Medicom Healthcare Ltd.,
Lynton House, 7-12 Tavistock square,
Kings Cross, London, WC1H 9LT, U.K.

Mfg. Lic. No.: MB/06/470

7-HPMC.3SDMD-19-GB-PV5